




# MS Microchip Syringe

## 0409999-10 MS Microchip Syringe 2x 12 mm

MS Schippers offers a highly ergonomic syringe design that facilitates the application of the microchip. The syringe has a locking system so the microchip is not applied or dislodged accidentally. It makes a "click" sound when the microchip has been applied reassuring the veterinarian that it has been implanted adequately. The size (2x 12 mm) of the transponder is small enough that it permits subcutaneous implantation and does not harm the animal. Implantation is done by using a sterilized hypodermic needle. The microchip upon insertion is impossible to lose, to be altered or to be removed. It is coated in a product called parylene so it does not move throughout the body of the animal after application, therefore it remains in the same place during the life of the animal. The microchip's memory contains a unique number. Two animals can never have the same code. In addition, Schippers produces the microchip FDX-B technology, according to the ISO standard 11784/11785. MS Microchip Syringes are ICAR-registered.

## Packing and Presentation

*This product is supplied in boxes that contain 10 units. The identification unit is supplied in a sterilized blister with 6 barcode labels, a needle and a syringe.*

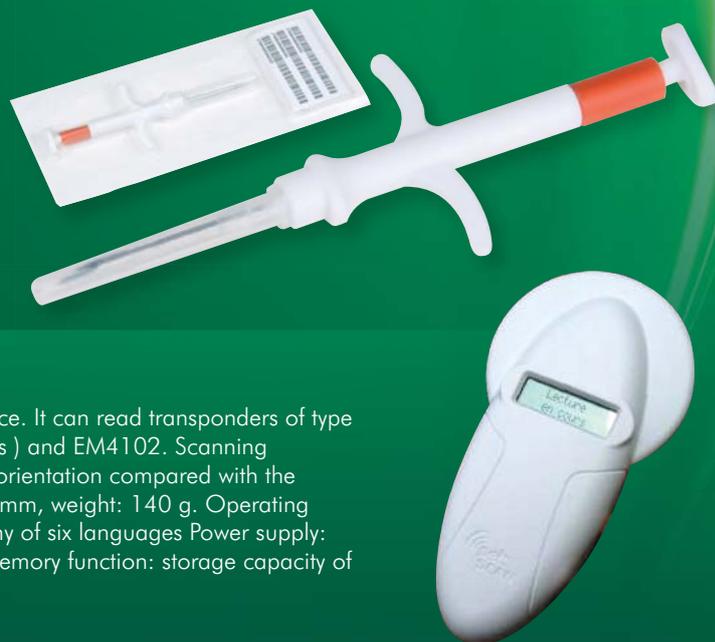

## 0402575 Petscan Reader

The transponder reader requires little space but has excellent performance. It can read transponders of type FDX B (ISO standard), type FDX A (Destron, AVID, AEG), HDX (Texas Tiris ) and EM4102. Scanning distance: from 5 to 8 cm. It depends on the type of transponder and its orientation compared with the scanner. Dimensions and weight: L: 155 mm, W: 80 mm, thickness: 33 mm, weight: 140 g. Operating temperature: - 0°C to + 40°C Display: 2 lines of 16 characters using any of six languages Power supply: 9V alkaline or lithium battery of type PP3. Do not use saline batteries. Memory function: storage capacity of 1000 identifiers transmissible to a PC or PDA via the RS 232 link.

**MS SCHIPPERS**  
Passion for farming

## Schippers Export BV

Rond Deel 12 - 5531 AH Bladel • The Netherlands  
Tel: +31(0)497-339774 • Fax: +31(0)497-339779  
export@msschippers.com • www.msschippers.com